

Aung Myin Kyaw (Max) weekly Research Progress Report

Date : 18/08/2018 to 25/08/2018

➤ Scope of the work:

- ✓ Printing of the data to csv file and understanding the values**

➤ **Research Progress in last week**


- This is week I was able to collect data into the csv file (the file is attached as 'test1.csv'). The code I used to collect data is attached in the following slide.
- However, the data collected does not tally with the data I downloaded from Thorlabs which is attached as 'testthorlabs.csv'. Thus, for the coming week I will be examining the representation of the values that are collected.

```
CCS100.startScan();
CCS100.getScanData(Data);
Console.WriteLine(Data.Length);
StreamWriter sw = new
StreamWriter(@"C:\Users\AungMyin\Desktop\School\FYP\Final-Year-
Project\test1.csv", true);
for (int i = 0; i < Data.Length; i++) {
    {
        sw.WriteLine(Data[i].ToString());
        //sw.Write("\n");
    }
}
sw.Write("\r\n");
sw.Flush();
sw.Close();
```

By using the function that I presented above, I am able to save the scan data by CS100 to excel file.

➤ Research Progress in last week

- In addition, from this experiment I am able to find out that minimum (321.896381278245) and maximum (746.289219694925) wavelength that can be measured using my wavelength. The wavelengths are the same when I compared with the maximum and minimum wavelength that measured using Thorlabs.

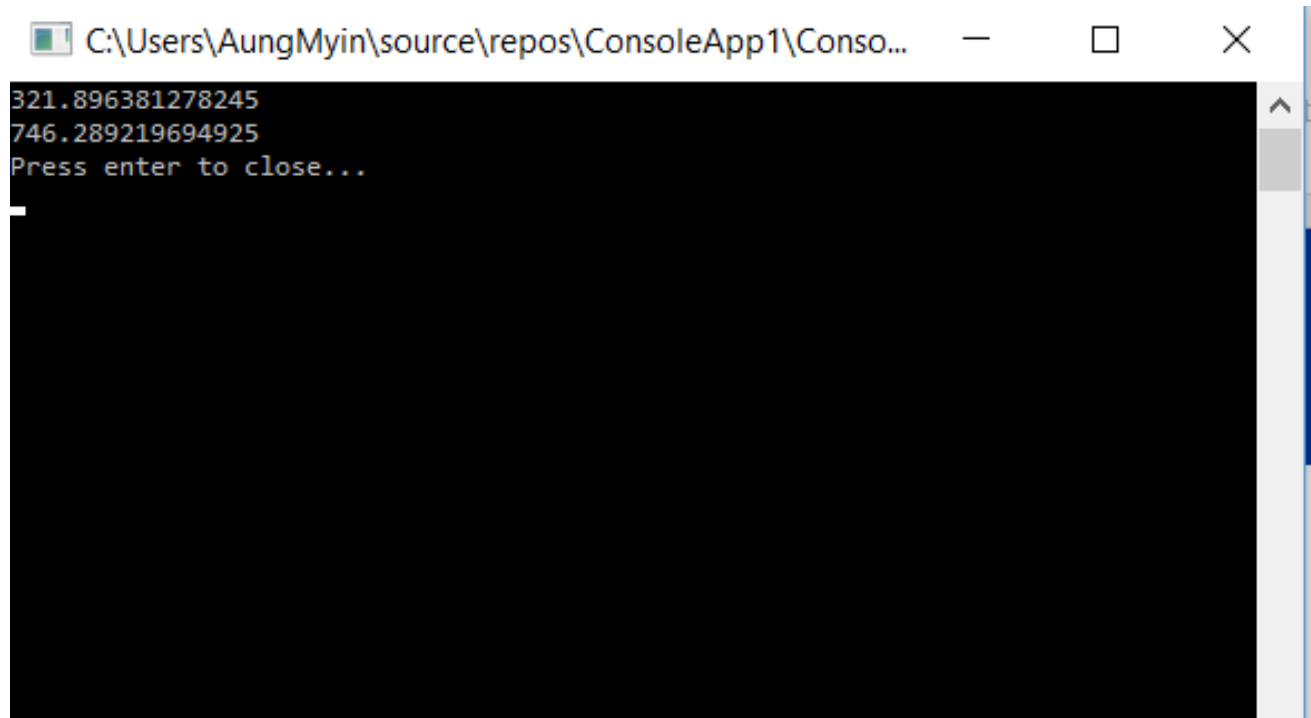

```
C:\Users\AungMyin\source\repos\ConsoleApp1\ConsoleApp1.exe
321.896381278245
746.289219694925
Press enter to close...
_
```

➤ **Research Plan for next week**

- ❖ Identifying the values that are being save to the csv file. Check whether any calibration is needed to collect the data correctly.
 - ❖ Once the values are identify, try to read those values and plot a graph in the program.
- .